

BOLD[®]

COMPRESSION SCREW



- **Functionality - Small Bone Fragment Fixation**
- **Flexibility - Comprehensive Product Portfolio**
- **Stability - Titanium Implant**

**LOWER
EXTREMITY
SOLUTIONS**

 **INTEGRA**[™]
Extremity Reconstruction

BOLD® COMPRES




Cannulated screw




Lot number and reference* marking for easier product traceability
* except size 10 mm

FOREFOOT SET

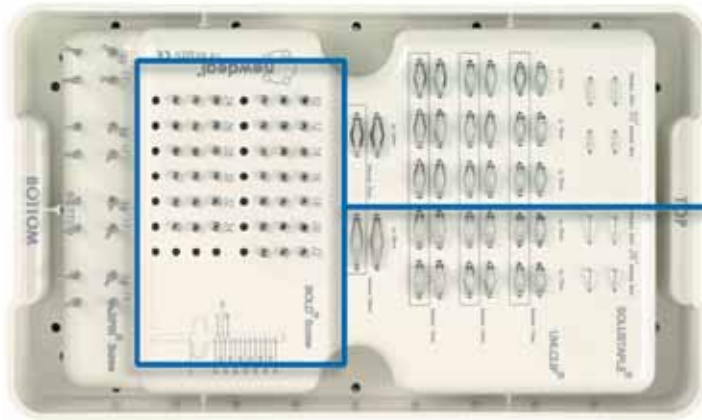

BOLD SCREWS

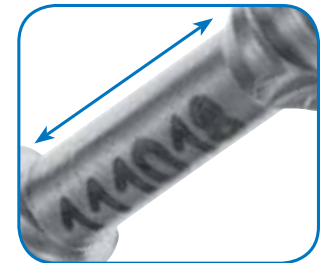

Variable lag part




BOLD® CLAMP

DRILL DIAM 2.2 MM


K-WIRE LENGTH 70 MM

«2 IN 1» DRILL

BOLD® SCREWDRIVER

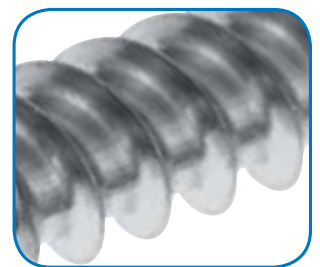

Dual threads to increase grip and create optimal compression

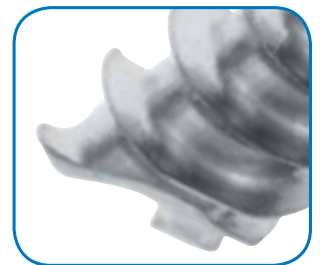

Self-tapping tip with 3 notches to ease head introduction


SION SCREW



Description

Titanium compression screw diameter 3.0 mm.

- Totally intra-osseous implant.
- Large range of size: 10 to 34 mm by 2 mm increment.
- K-Wire guided drilling and screwing.
- Self-tapping.
- True & Controlled compression.
- In order to optimize compression, the head has to be totally embedded.
- Material: Titanium alloy
TIAI6V4, ISO 5832-3 ASTM F136.

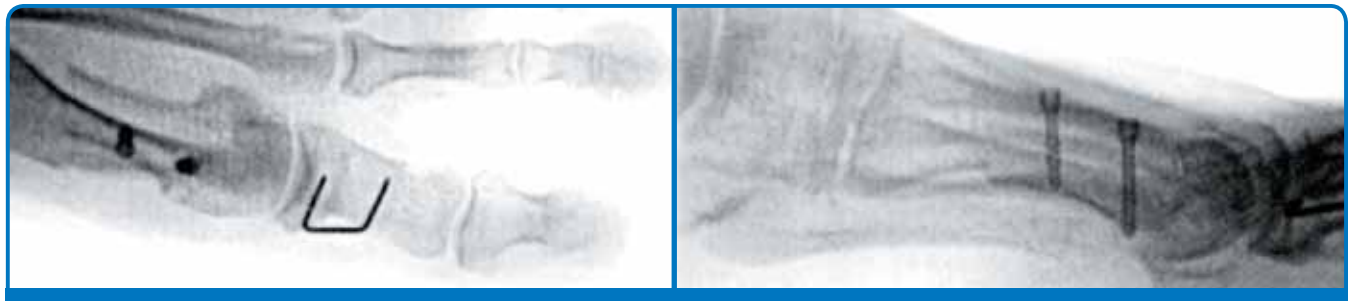

INDICATIONS*

The Bold screw is indicated for fixation of bone fractures or for bone reconstruction.

- Fixation of small bone fragments.
- Arthrodesis in foot surgery.
- Mono or bi-cortical osteotomies in the foot.
- Distal or proximal metatarsal osteotomies.
- Fixation of osteotomies for Hallux Valgus treatment (such as Scarf, Chevron, etc.).

OTHER POSSIBLE INDICATIONS

- Arthrodesis in hand surgery.
- Mono or bi-cortical osteotomies in the hand.
- Distal or proximal metacarpal osteotomies.

* See package insert for full prescribing information.

BOLD®

COMPRESSION SCREW

BOLD Diam. 3.0mm Compression Screw

Catalog Number	Length	Catalog Number	Length
111 010ND	10 mm	111 022ND	22 mm
111 012ND	12 mm	111 024ND	24 mm
111 014ND	14 mm	111 026ND	26 mm
111 016ND	16 mm	111 028ND	28 mm
111 018ND	18 mm	111 030ND	30 mm
111 020ND	20 mm	111 032ND	32 mm
		111 034ND	34 mm

Instrumentation

Catalog Number	Description
115 070ND	K-Wire length 70 mm
119 103ND	Cannulated measurer screwdriver
119 105ND	Bold clamp
119 004ND	Cannulated drill diam. 2.2 mm
119 005ND	«2 in 1» drill length 18 mm



- The products are manufactured and referenced within the frame of the standards in force.
- Implantation procedures are described in the surgical technique.
- WARNING: Federal law (USA) restricts this device to sale by or on the order of a physician.

THE FOREFOOT SET ALSO INCLUDES THE FOLLOWING PRODUCTS:




UNI-CLIP®



SOLUSTAPLE®



BOLD®

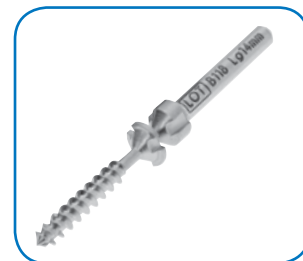

SPIN®

 **INTEGRA™**
Extremity Reconstruction

Integra LifeSciences Corporation
311 Enterprise Drive • Plainsboro, NJ 08536
877-444-1122 (USA) • 609-936-5400 (Outside USA) • 866-800-7742 (Fax)
www.Integra-LS.com

